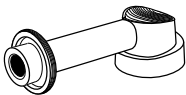
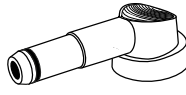


Maintenance instructions

<p>Cleaning preparation</p>	<p>Remove soiling from the electric motor with a disposable cloth/paper towel. Unscrew the motor cap and remove the cable. For the Electronic Motor 21 additionally remove the handpiece carrier.</p> <p>Electronic Motor 21</p> <p>Electronic Motor 31ESS</p>
<p>Inspection and maintenance</p>  <p>Art. No. 1974</p>  <p>Art. No. 1942</p>	<p>Perform a visual inspection to check for damage, corrosion, wear and tear. Connect Electronic Motor 21 with the spray nozzle (1974) by screwing it on and spray motor for about 3 seconds. For the Electronic Motor 31ESS the nozzle (1942) is plugged in from the front onto the handpiece carrier. Briefly wipe with a moist cloth (see instructions on spray can). After spraying, screw the handpiece carrier and the motor cap with cable back onto the electronic motor.</p> <p>Electronic Motor 21</p> <p>Electronic Motor 31ESS</p>

Service Centers

Switzerland
 Nouvag AG • St.Gallerstr. 23-25 • CH-9403 Goldach
 Phone +41 (0)71 846 66 00 • Fax +41 (0)71 845 35 36
 info@nouvag.com • www.nouvag.com

Germany
 Nouvag GmbH • Schulthaißstrasse 15 • D-78462 Konstanz
 Phone +49 (0)7531 1290-0 • Fax +49 (0)7531 1290-12
 info-de@nouvag.com • www.nouvag.com